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Modified Headgear

Casques ajustables

Anclaje universal

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Instructions for use

Mode d'emploi

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The Hickham Modified Headgear

Dear Customer

Thank you for choosing this high-quality product from Dentaureum.

In order to achieve the best results with this product, it is important to carefully study and follow these instructions for use.

The written instructions cannot cover all possible situations that can occur during use. For this reason, our hotline is available to answer any other questions or concerns that may arise.

Due to constant developments, we recommend, even when using the same products regularly, that you study the enclosed current instructions for use or refer to our website at www.dentaureum.de.

1. Manufacturing firm

Dentaureum GmbH & Co. KG · Turnstraße 31 · 75228 Ispringen · Germany

2. With reference to our quality standards

Dentaureum ensures a faultless quality of the products manufactured by us. These recommendations are based upon our own experiences. The user himself is responsible for the processing of the products. Responsibility for failures cannot be taken, as we have no influence on the processing on site.

3. General information

This headgear provides the orthodontist with simplified fitting of headgears to patients since only one basic model is required for all head-sizes whilst the pull direction can be individually selected. The required forces are generated by rubber rings of different strengths.

Bibliography

Hickham, J. H. and Miethke, R. R.:
Data on Jaw Orthopedics and Orthodontics, Vol. 4, 1979,
The Quintessenz, Vol. 10, 1979.

A. Horizontal Pull Headgear, adjustable

Range of applications

Intrusion of individual teeth: For this, the J-hooks (REF 744-713-00) are attached mesially over the archwires directly against the canines or premolars.

Intrusion of groups of teeth: Control of the growth of the upper and lower jaws.


Force direction – see illustration.

Fitting

To fit the patient correctly, the plastic guides are adjusted with both hands:

- a) vertical:
The middle of the ears and the centers of the plastic guides must coincide.
- b) sagittal:
The plastic guides must be positioned as close as possible in front of the ear so that the forces are applied to the bony support (os temporalis).

The plastic guide should be tilted slightly forwards so that the distance between the perforations of the guide and the corner of the mouth is constant. This will lead to the headgear/plastic guide rivets lying perpendicular over one another.

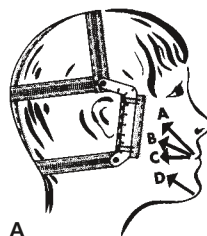
A

REF 744-712-00

1 piece

Horizontal pull
headgear including
cheek pads


Correct positioning of the plastic guides provides for the correct fit of the headgear to the patient and its being held firmly in place.

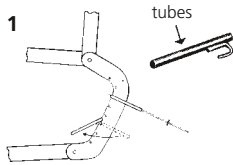Attachment of the tubes (fig. 1)

The hook of the tube is placed in one of the open holes on the inner side of the plastic guide and angled downwards. Do not force it through the slot.

Lift the tube and without using force turn it forwards by 180°. Reverse this procedure if the tube has to be replaced or its position changed.

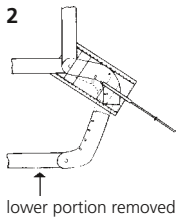Attachment of the J-hooks (REF 744-713-00)

Push the J-hook into the tube and attach one of the previously selected rubber rings to the double hooks and to the far end of the tube.

Extraoral latex elastics, standard size = 9.5 mm/ $\frac{3}{8}$ " (REF 772-205-01)

Variations (fig. 2)

In the later stages of orthodontic treatment, only one pair of J-hooks in the high pull position is needed. In this case, it is usually an advantage to remove the lower part (the neck portion) of the headgear. For this, the plastic guide is cut in half with side cutters and the vertical strap is also cut through behind the head just below the horizontal strap. The cheek pads should then be slid over the upper part of the guide and stapled in position. The headgear will now look like a conventional high-pull appliance.

B. High Pull Mode, adjustable

Fitting


Fitting of this headgear is very simple since there is only one ready-for-use size and no further stapling is needed. The headgear is placed on the patient and both plastic guides (together with cheek pads and tubes) are drawn upwards or downwards as required. A stud-anchor mechanism simplifies fixation in the desired position. Five fixation positions are possible.

The exact position is found by the patient holding the J-hooks, which are anchored to the inner bow, in the chosen pull direction.

The two plastic strips are adjusted to the same length with the distance between the double hooks (of the J-hooks) and the start of the tubes being about 1.5 – 2.0 cm.

Changing the cheek pads

In order to change the cheek pads, the J-hooks are removed, the old cheek pads taken off and new ones slid over the tubes – see illustration.

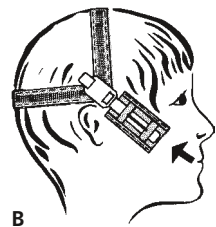B

REF 744-717-00


1 piece

High pull headgear including cheek pads

Accessories

J-hooks, complete with tubes and rubber rings. For distal movement of canines/premolars.

For Universal pull headgear (REF 744-712-00) and high pull headgear (REF 744-717-00).

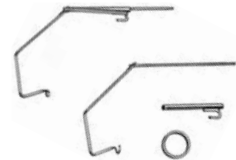REF 744-713-00

2 pieces

Extraoral latex elastics

A range of elastics for all extraoral techniques.

The stated elastic pull comes into force when the elastic is stretched to three times its diameter.


Standard size 9.5 mm/ 3/8"

Elastic pull	Inner diameter		REF 1000 pieces each
	mm	inch	
Light pull 227 g = 2.3 N	6.4	1/4"	772-203-01
	8.0	5/16"	772-204-01
	9.5	3/8"	772-205-01
	12.7	1/2"	772-206-01
Heavy pull 455 g = 4.5 N	8.0	5/16"	772-208-01
	12.7	1/2"	772-210-01

Headgear Adjustment

Position the headgear on the patient's head. To ensure a good fit, push the loose ends of the straps through the two-way securing element (fig. A) or locking element (fig. B) as applicable.

Using the stapler, secure the straps with 3 staples (fig. C) and cut off ends.

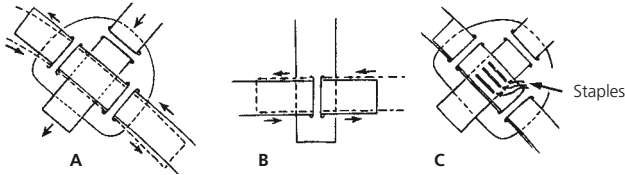
Chin pad

For lining the chin cap
Recommended for use in high temperatures, high humidity, or if the patient suffers from allergies.

REF 744-616-00

In order to avoid injury, the open ends of the staples must point away from the head.

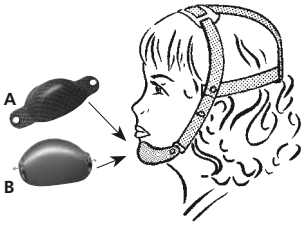"High Pull" headgear Designed for chin-cap treatment

A = with chin cap (elastic, porous, air permeable)


REF 744-705-00

B = with chin cap (rigid with small external hooks)

REF 744-709-00

Headgear color = blue

One type is suitable for all head and chin sizes.

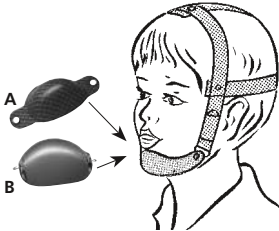"Vertical Pull" headgear Designed for chin-cap treatment

A = with chin-cap (elastic, porous, air permeable)

REF 744-706-00

B = with chin cap (rigid with small external hooks)

REF 744-710-00

Headgear color = blue

One type is suitable for all head and chin sizes.

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